

Escape-Carry Chair® ST / XS / ST-PLUS / CF

Evacuation ascending the stairs, while in a sitting position.

What's an Escape-Carry Chair®?

An Escape-Carry Chair® is the solution to evacuate people to safety in a seated position up and/or down a flight of stairs.

A simple procedure can make the Escape-Chair® fully functional within a few seconds.

How does an Escape-Carry Chair® work?

The unique glide/brake system makes the braking of the Escape-Carry Chair® possible on all different kinds of stairs, from wood to concrete. This braking system is optimized by the weight of the evacuee, in combination with the pressure of the operator pushing down on the handle in the direction of the stairs.

Ascending evacuation should always be operated by at least 2 people and ensure to make use of the support braces.


Advantages of an Escape-Carry Chair®:

- Used for ascending and/or descending evacuation.
- Ideal for quick evacuation.
- Descending evacuation can be operated by one person.
- Speed of descent to be determined by the user (extra brake is unnecessary).
- Ready for use within seconds.
- Easily maneuvering into the evacuation flow.
- Forms minimal blockage on the stairs.
- Equipped with carrying handles.
- Fits into your emergency plan.
- Safe closed gliding belt cassette.
- CE, TÜV & ISO 9001 certified.

Our advice:

Make sure you have an Escape-Carry Chair® available in your organization so that you're prepared for any emergency, for up and/or down evacuation.

"Failure to prepare for an emergency has its consequences."

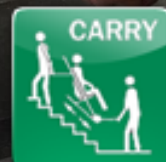

ST Model

XS Model

ST-PLUS Model

CF Model



Specifications	Escape-Carry Chair® ST	Escape-Carry Chair® XS	Escape-Carry Chair® ST-PLUS	Escape-Carry Chair® CF
Article number:	10110	10150	10120	10130
Dimensions folded up (height x width x depth):	102 x 51 x 21 cm	105 x 40 x 27 cm	120 x 51 x 21 cm	120 x 57 x 21 cm
Maximum payload:	185 kg	185 kg	185 kg	185 kg
Maximum stair gradient:	40°	40°	40°	40°
Weight:	±13 kg	± 12 kg	± 14 kg	± 16 kg
Ergonomische backrest:	single	single	double	double
Demonstration and advise provided:	✓	✓	✓	✓
Ready for use within seconds:	✓	✓	✓	✓
Speed of descent to be determined by the user:	✓	✓	✓	✓
Safe closed gliding belt cassette:	✓	✓	✓	✓
Reflective strips:	✓	✓	✓	✓
CE, TÜV & ISO 9001 certified:	✓	✓	✓	✓
Delivery incl. 4-language manual, instruction DVD, identity sign & wall hooks:	✓	✓	✓	✓
Footrest:	optional	-	optional	✓
Armrests:	optional	-	optional	✓
Padded (upholstered) headrest:	optional	✓	✓	✓
Comfort seat:	-	-	✓	✓
Two height adjustable operating handle:	-	-	✓	✓
Extra slim design (e.g. for airports):	-	✓	-	-
Warranty on non-wearing parts:	3 years	3 years	3 years	3 years



Interim changes to technical specifications and printing errors reserved.
Produced by Escape Mobility Production & Innovation: TÜV approved organization.

