



Evacuation directly ... from a bed

What`s an Escape-Sheet®?

People who aren't able to leave their bed independently need to be evacuated using evacuation equipment. One of these aids is the Escape-Sheet®. A mattress with an Escape-Sheet® underneath can be transformed very quickly into an evacuation mattress.

How does an Escape-Sheet® work?

In case of emergency the strong safety belts are taken from underneath the mattress. The belts are strapped over the patient and bedding so the patient is safely secured to the mattress. The Escape-Sheet® can then be easily pulled to the floor and taken to a place of safety.


Advantages Escape-Sheet® HFR:

- HFR design: heat and flame-resistant (EN ISO 11612).
- Static-Control™ (anti-static polyester).
- Hydro-tech Nano technology (water-repellent)
- Straps and handles highly visible in bright green colour.
- Extra-long towing straps to help operators with improved handling practices.
- Plastic buckle pockets increase the life expectancy of the product (allocated on the side of the sheet).
- Handgrips to assist with sideways manoeuvrability.
- Adjusting corner straps (suitable for each mattress height).
- TÜV & CE certified.

Our advice:

Make sure every bed in your Organization is equipped with an Escape-Sheet®, so you're always prepared for optimal and rapid evacuation.



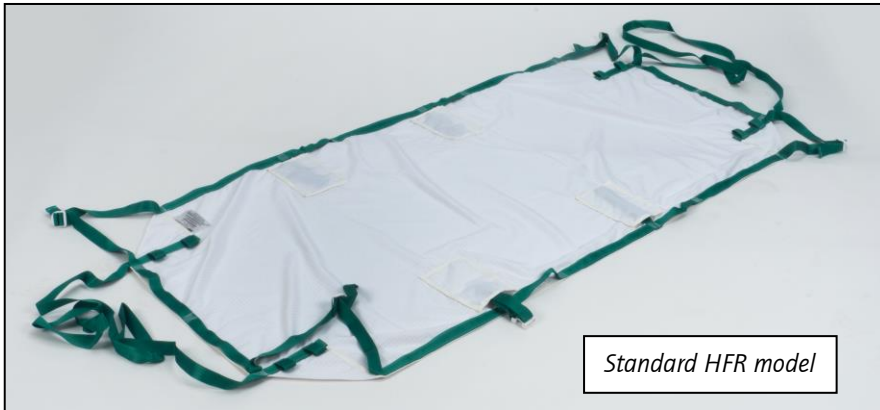

Standard HFR model



Option: Straps in 3 different colors



Extra-long towing straps



Corner sheet adjusters




Buckles in pockets



Washing instructions

Specifications	Escape-Sheet® HFR Standard	CE
Article number:	21241	
Dimensions (length x width):	± 200 x 90 cm	
Material:	54/45/1% PPAN-fr/lyocell/Static-Control™	
Maximum payload:	200 kg	
Weight:	± 1100 grams	
Demonstration and advise provided:	✓	
Ready within seconds:	✓	
Immediate evacuation directly from the bed:	✓	
Heat and flame-resistant (EN ISO 11612):	✓	
Industrial wash 75° C:	✓	
Hydro-tech nano technology (water-repellent):	✓	
Adjustable corner straps:	✓	
For all mattress thicknesses:	< 23 cm	
Buckles in pockets:	✓	
Towing straps available in 3 different colors (intuitive handling):	Optional	
Warranty:	2 years	

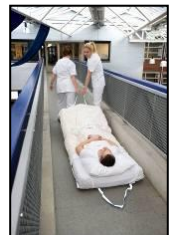

Interim changes to technical specifications and printing errors reserved.
 Produced by Escape Mobility Production & Innovation: TÜV approved organization.
 Health guaranteed through Öko-tex Standard 100

